



HeartSine®

samaritan™ PAD

Automated External Defibrillator

Rugged, Portable, Easy-to-Use

Lifesaving Technology for Utility Companies.



Sudden Cardiac Arrest (SCA) is a leading cause of death globally and response time is critical for survival. The HeartSine's samaritan PAD is durable enough to go with you anywhere, anytime and will operate under most any environmental condition.




Unmatched Durability

Portable. The samaritan PAD is lighter (2.4 lbs.) and smaller than other defibrillators.

Unmatched durability. The samaritan PAD carries the industry's highest durability rating – an **IP56**. With this rating, it can be used in virtually any environmental condition and will withstand **dust** and **water jetting**. In addition, it will withstand **Shock** (MIL STD 810K Method 516.5) and **Vibration** (Method 514.5). It also carries a **7-year unit warranty**.

Advanced technology. SCOPE™ Biphasic technology automatically optimizes energy output for each person, and provides significantly better performance in removing ventricular fibrillation (VF) by the third shock.**

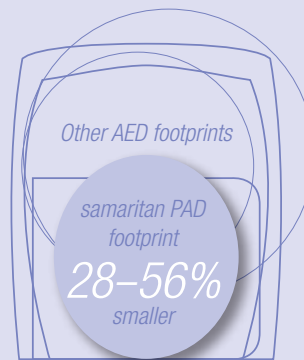

2.4 lbs. light.

Easy-to-follow visual and verbal guides.

User-friendly. The samaritan PAD features easy-to-understand visual and oral prompts that guide a user through the process.

Two-button operation. Only two buttons, ON and SHOCK, are required, providing straightforward operation.

Always ready. A System Status Ready Indicator flashes to show that the complete system is operational and ready for use. Device automatically runs self check each week.



Visual cues prompt pad placement



Stand clear of the patient



It is safe to touch the patient

Real economy for the real world.

Two parts, one expiration date. Pad-Pak™ cartridge combines battery and electrode pads, with one expiration date to monitor.

Low cost of ownership. Cartridge has a shelf-life of 3.5 years from date of manufacture, offering significant savings over other defibrillators that require separate battery and pad units.



Pad-Pak™ and Pediatric-Pak™ with pre-attached electrodes



Technical Overview

Physical

Size: 8.0 in x 7.25 in x 1.9 in (20cm x 18.4cm x 4.8cm)
Weight: 2.4 lbs (1,1 kg) including Pad-Pak™ Battery

Defibrillator

Waveform: SCOPE® (self-compensating output pulse envelope) Biphasic waveform. Optimized biphasic waveform compensates energy, slope and envelope for patient impedance

Energy Selection

Pre-configured factory settings for escalating energy are AHA/ERC 2005
Adult
1. Shock 150J 2. Shock 150J 3. Shock 200J
Paediatric
1. Shock 50J 2. Shock 50J 3. Shock 50J

Patient Analysis System

Method: Evaluates patient's ECG, signal quality, electrode contact integrity and patient impedance to determine if defibrillation is required
Sensitivity/Specificity: Meets AAMI DF80:2003 requirements

Environmental

Operating Temperature: 0°C to 50°C (+32°F to +122°F)
Storage Temperature: -10°C to 50°C (+14°F to +122°F)
Relative Humidity: 5% to 95% (non-condensing)
Water Resistance: IEC 60529/EN 60529 IP56
Altitude: 0 to 15,000 feet (0 – 4,575 meters)
Shock: MIL STD 810F Method 516.5, Procedure I (40G's)
Vibration: MIL STD 810F Method 514.5+
Category 4 Truck Transportation – US Highways
Category 7 Aircraft – Jet 737 & General Aviation (Exposure)
EMC: EN 60601-1-2, Second Edition: 2001
Radiated Emissions: CISPR11: 1997 +A1:1999+A2:2001 Group 1 Class B
Electrostatic Discharge Immunity: EN61000-4-2:1995 (8KV)+A1:1998+A2:2000
RF Immunity: EN61000-4-3:1996, 80 MHz – 2.5 GHz, (10V/m)+A1:1998+A2:2000
Magnetic Field Immunity: EN61000-4-8:1993 (3 A/m)+A1:2000
Aircraft: RTCA / DO – 160D: 1997, Section 21 (Category M)
Falling height: 1 meter

Event Documentation

Type: Internal memory
Memory Capacity: 45 minutes of ECG (full disclosure) and event/incident recording
Playback Capabilities: Custom USB cable directly connected to PC and Saver™ EVO Windows-based data review software

Materials Used

PAD SAM300P: ABS, Santoprene. Printed circuit board with electronic components.
PAD Cartridge: Battery: Lithium Manganese Dioxide
Housing: ABS – Electrodes: Hydrogel, Silver, Aluminium and Polyester

Combined Electrode and Battery Cartridge

Pad-Pak™
Shelf Life: Check Expiration Date (typically 3.5 years from manufacture)
Weight: 0.44 lbs (0.2kg)
Size: 3.93 in x 5.24 in x .94 in (10cm x 13.3cm x 2.4cm)
Battery Type: Lithium Manganese Dioxide (LiMnO₂)
18V, 0.8 Amp Hrs
>30 shocks at 200J or 6 hours of continuous monitoring.
Capacity:

Adult Electrodes (Pad-Pak-01)

Standby Life: Typically 3.5 years from manufacture date
Electrodes: Samaritan® disposable defibrillation pads are supplied as standard with each device
Anterior-lateral
Placement: 100cm²
Active Gel Area: 3.5 ft (1m)
Cable Length:

Children Electrodes (Pad-Pak-02)

Standby Life: Typically 3.5 years from manufacture date
Electrodes: Samaritan® disposable defibrillation pads are supplied as standard with each device
Anterior-Posterior
Placement: 100cm²
Active Gel Area: 3.5 ft (1m)cable
Cable Length:



*Self Compensating Output Pulse Envelope technology automatically compensates energy, slope and pulse envelope for the patient.

**Efficacy of Distinct Energy Delivery Protocols Comparing Two Biphasic Defibrillators for Cardiac Arrest, Walsh, McClelland, Owens, Anderson, Turner, Adgey; *The American Journal of Cardiology*, Vol. 94, Aug.1, 2004

Corporate Headquarters:

HeartSine Technologies, Inc.
105 Terry Drive,
Newtown, PA. 18940
Tel: (215) 860 8100
Fax: (215) 860 8192
Web: heartsine.com

Manufacturing Facilities

HeartSine Technologies, Ltd
Canberra House, 203 Airport Road West
Belfast, Northern Ireland BT3 9ED
Tel: +44 (0) 28 9093 9400
Fax: +44 (0) 28 9093 9401

CAUTION: U.S. Federal law restricts this device to sale by or on the order of a licensed practitioner.



The products described in this brochure all meet the applicable European Medical Directive requirements.